---

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Hongyi Wu(吴鸿毅)

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---

## Contents:

---

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The screenshot shows the Digitizer 2.0 Home page in Mozilla Firefox. The title bar reads "Digitizer 2.0 - Home — Mozilla Firefox (于 wuhongyi-u20)". The main content area displays board information for DT2745B [PID: 25480] and a list of errors flags.

**Board information:**

Model name:	DT2745B
Firmware version:	2022092904
Firmware type:	DPP_PSD
Bitstream version:	0.1.27
PID:	25480
License:	Not Licensed
Trial time left:	02:16
BIOS:	2.1
System:	2022090600

**Hardware monitor:**

Temp Sens Air In:	31.2 °C
Temp Sens Air Out:	44.4 °C
Temp Sens Core:	69.2 °C
Temp Sens First ADC:	40.3 °C
Temp Sens Last ADC:	54.1 °C
Temp Sens DC/DC:	43.9 °C
V In Sens DC/DC:	5.039 V
V Out Sens DC/DC:	0.725 V
I Out Sens DC/DC:	7.484 A
Speed Sens Fan 1:	3330 rpm
Speed Sens Fan 2:	3300 rpm

**Errors flags:**

Power Status	Green
Board Init	Green
SI5341 PLL Lock	Green
SI5395 PLL Lock	Green
LMK04832 PLL Lock	Green
JESD204B Lock	Green
DDR4 PL Bank0 Calibration	Green
DDR4 PL Bank1 Calibration	Green
DDR4 PS Calibration	Green
FPGA Configuration	Green
BIC Check	Green
ADC Temp	Green
Air Outlet Temp	Green
FPGA Temp	Green
DC/DC Temp	Green
Clock In	Green
ADC Shutdown	Green

**Left sidebar:**

- Menu (Dashboard selected)
- Board ID card
- Network settings
- User management
- Firmware
- Clock settings
- License
- Mon, 17 Dec 2103  
05:21:34 GMT
- Version 1.1.7  
©2021-2024 CAEN SpA

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The screenshot shows the 'Firmware' section of the GDDAQ web interface. On the left is a sidebar with various menu items: Dashboard, Board ID card, Network settings, User management, Time/Date settings, Parameters (BETA), Firmware (selected), Clock settings, License, and footer information (2024/2/6, GMT+1 14:29:01, Version 1.2.0, ©2021-2024 CAEN SpA). The main content area is titled 'Firmware' and contains a table of available firmware versions:

Version	Type	Status	Action Buttons
2023112703	ETH IG	Current	More
DPP-PHA			More
2023091902	ETH IG		More
DPP-ZLE			More
2023091901	ETH IG		More
DPP-PSD			More
2023091900	ETH IG		More
Scope			More
2022052302	ETH IG	Factory	More
Factory firmware Scope			More

A green progress bar at the top indicates 'Storage Space Usage: 11%'. A blue button labeled 'Upload new firmware...' is located above the table.

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- **ADC\_DATA**
  - Data from the ADC (normal operating mode)
- **ADC\_TEST\_TOGGLE**
  - Toggle between 0x5555 and 0xAAAA (test mode)
- **ADC\_TEST\_RAMP**
  - 16-bit ramp pattern (test mode)
- **ADC\_TEST\_SIN**
  - 8-point sine wave test pattern
- **ADC\_TEST\_PRBS**
  - 16-bit PRBS generated by a 23-bit PRBS pattern generator (test mode)
- **Ramp**
  - Data from a ramp generator. It is actually a 16-bit field, where the 6 most significant bits identify the channel and the 10 less significant bits are the samples of a ramp from 0x000 up to 0x3FF (i.e. 0 to 1023). It is so a 10-bit ramp with offset given by “channel\*1024”. For channel 0, it is a counter from 0 to 1023; for channel 1, it is a counter from 1024 to 2047, and so on
- **IPE**
  - Not implemented
- **SquareWave**
  - Internally generated programmable square wave

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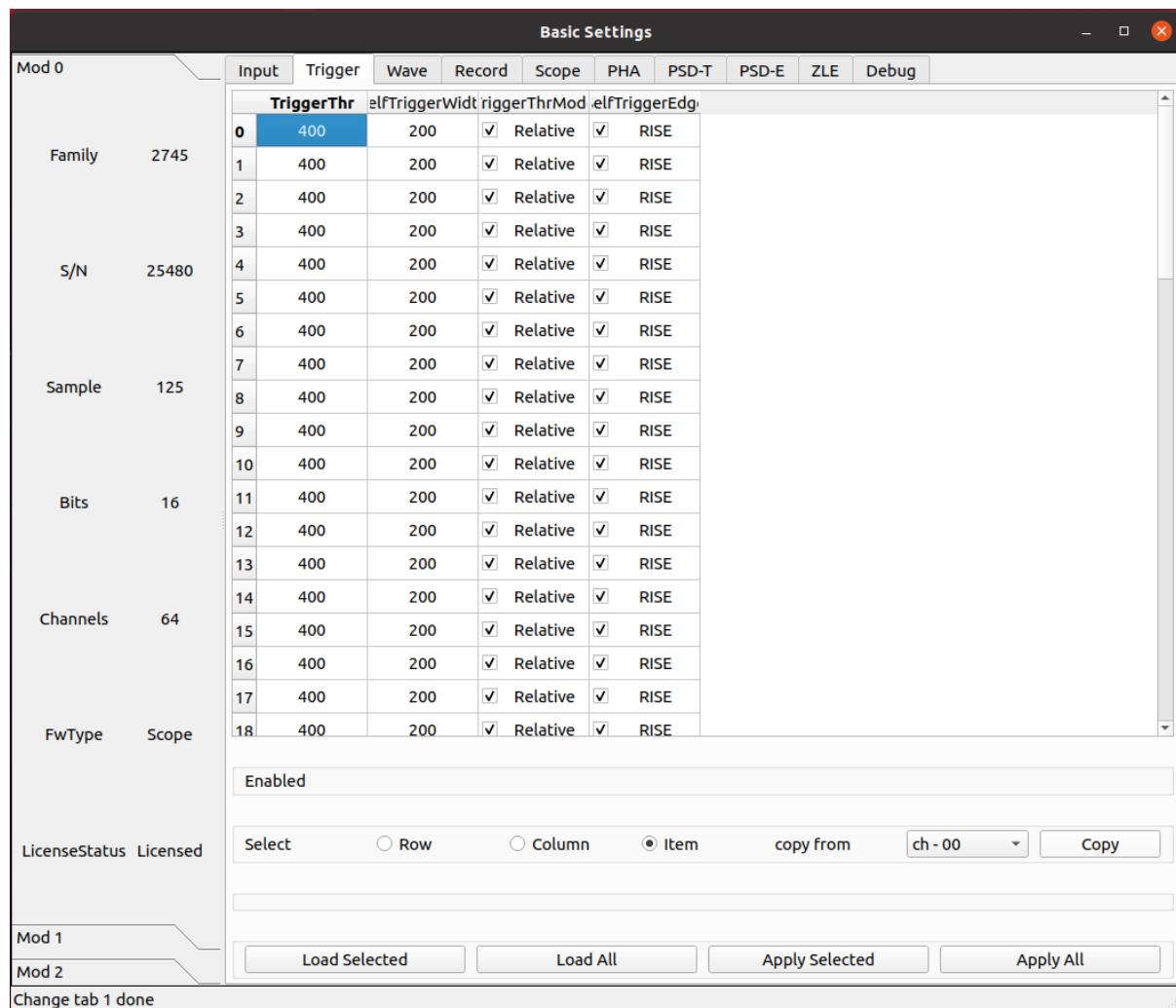
- **Relative**

- The threshold is relative to the baseline and automatically follows it when the DCOffset parameter is changed.

- **Absolute**

- The threshold is absolute, referred to the ADC input range. It does not follow the DCOffset setting.

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
- **RISE**

- The trigger occurs on the rising crossing of the threshold and the OverThr is true when the signal is above the threshold

- **FALL**

- The trigger occurs on the falling crossing of the threshold and the OverThr is true when the signal is below the threshold

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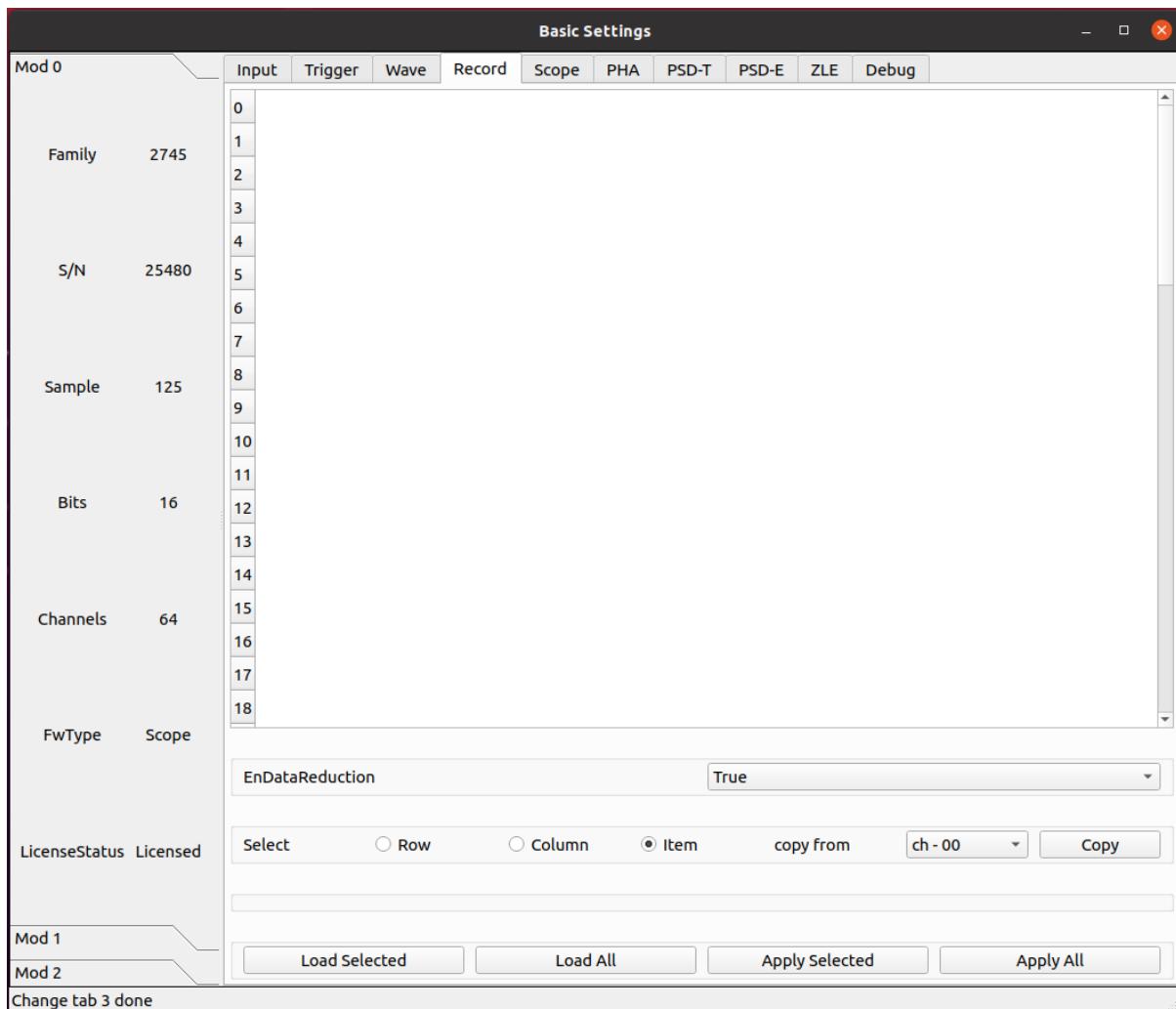
参数 EnDataReduction

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### 7.1.4 scope
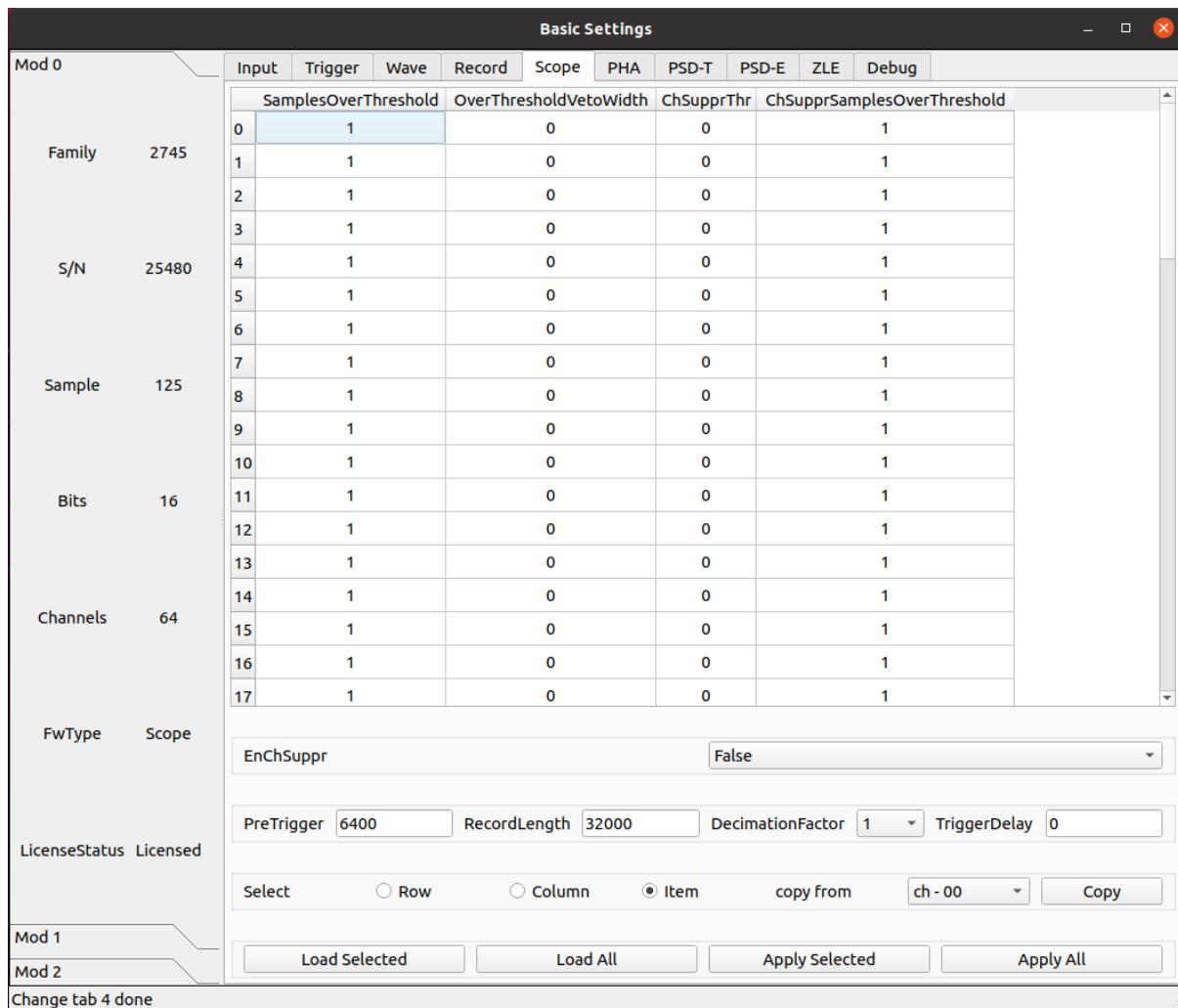
#### 参数 SamplesOverThreshold

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#### 参数 OverThresholdVetoWidth

Veto width to discard triggers when crossing the threshold in the opposite direction to the trigger one

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### 参数 EnChSupper

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- **True**
  - Channel zero-suppression is enabled
- **False**
  - Channel zero-suppression is disabled

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- **1**
  - Decimation disabled (default)
- **2**
  - Sampling Frequency / 2
- **4**
  - Sampling Frequency / 4
- **8**
  - Sampling Frequency / 8
- **16**
  - Sampling Frequency / 16
- **32**
  - Sampling Frequency / 32
- **64**
  - Sampling Frequency / 64
- **128**
  - Sampling Frequency / 128


- 256
  - Sampling Frequency / 256
- 512
  - Sampling Frequency / 512
- 1024
  - Sampling Frequency / 1024

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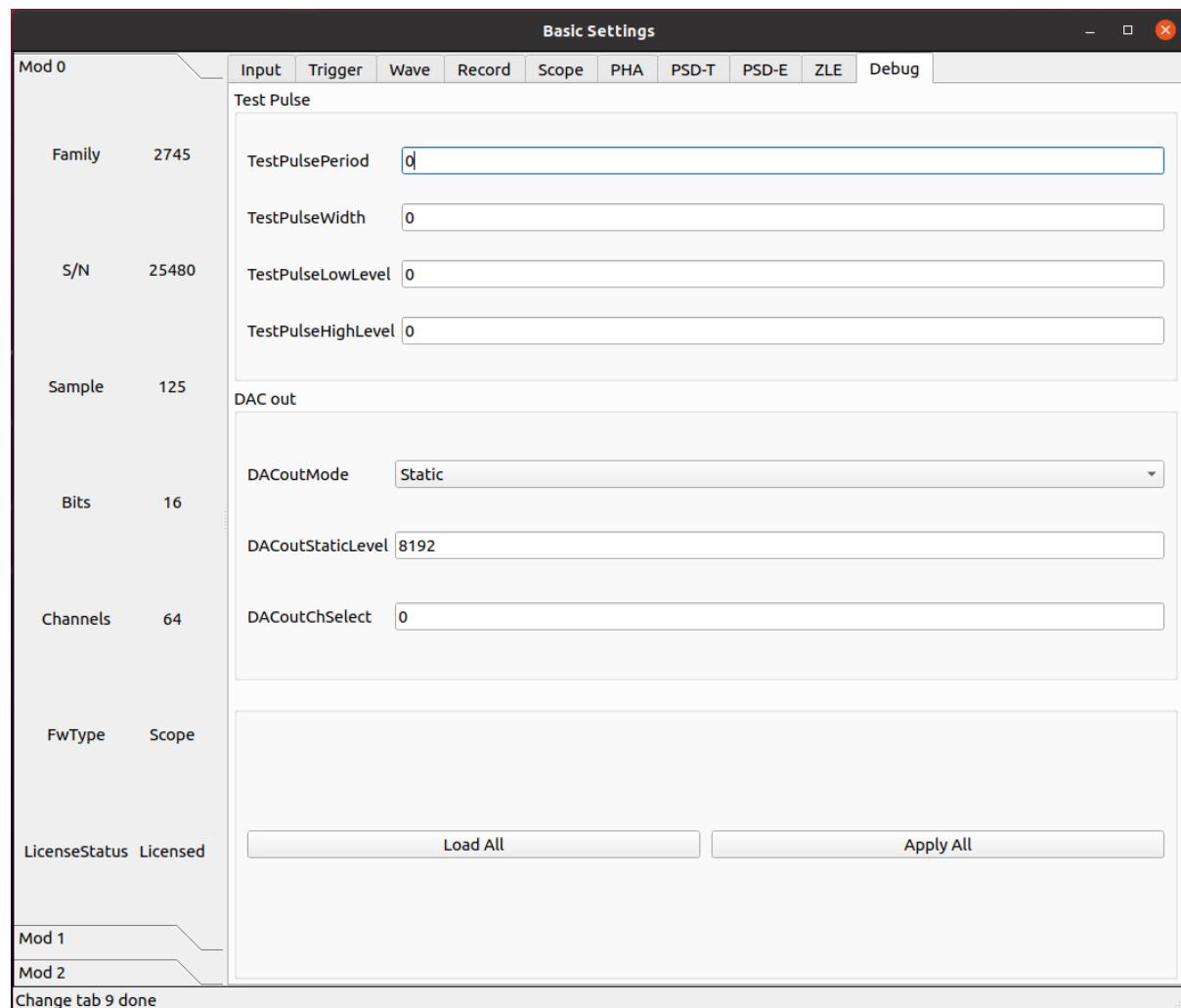- **Static**
  - DAC output stays at a fixed level, given by the DACoutStaticLevel parameter
- **Ramp**
  - The DAC output is driven by a 14-bit counter
- **Sin5MHz**
  - The DAC output is a sine wave at 5 MHz with fixed amplitude
- **Square**
  - Square wave with period and width set by TestPulsePeriod and TestPulseWidth and amplitude between TestPulseLowLevel and TestPulseHighLevel.
- **IPE**
  - Not implemented
- **ChInput**
  - The DAC reproduces the input signal received by one input channel, selected by the DACoutChSelect parameter
- **MemOccupancy**
  - Level of the memory occupancy (not yet implemented)
- **ChSum**
  - The DAC reproduces the “analog” sum of all the digitizer inputs (not yet implemented)
- **OverThrSum**
  - The DAC output is proportional to the number of channels that are currently above the threshold

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- **Disabled**
  - No trigger source enabled for the waveform
- **Ch64Trigger**
  - One (or more) channel self-trigger can generate a trigger for a waveform
- **ChSelfTrigger**
  - Channel self-trigger can generate a trigger for a waveform
- **SwTrg**
  - Software Trigger can generate a trigger for a waveform
- **ADCOverSaturation**
  - ADC Oversaturation can generate a trigger for a waveform
- **ADCUnderSaturation**
  - ADC Undersaturation can generate a trigger for a waveform
- **ExternalInhibit**
  - Inhibit can generate a trigger for a waveform
- **TRGIN**
  - External TRGIN can generate a trigger for a waveform
- **GlobalTriggerSource**
  - Acquisition Trigger Source (the same of the Scope mode) can generate a trigger for a waveform
- **LVDS**
  - A signal on the LVDS connectors can generate a trigger for a waveform
- **ITLA**
  - Internal Trigger Logic A can generate a trigger for a waveform
- **ITLB**
  - Internal Trigger Logic B can generate a trigger for a waveform

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- **Disabled**
  - No trigger source enabled for the T-E event
- **Ch64Trigger**
  - One (or more) channel self-trigger can generate a trigger for a T-E event
- **ChSelfTrigger**
  - Channel self-trigger can generate a trigger for a T-E event
- **SwTrg**
  - Software Trigger can generate a trigger for a T-E event

- **TRGIN**
  - External TRGIN can generate a trigger for a T-E event
- **GlobalTriggerSource**
  - Acquisition Trigger Source (the same of the Scope mode) can generate a trigger for a T-E event
- **LVDS**
  - A signal on the LVDS connectors can generate a trigger for a T-E event
- **ITLA**
  - Internal Trigger Logic A can generate a trigger for a T-E event
- **ITLB**
  - Internal Trigger Logic B can generate a trigger for a T-E event

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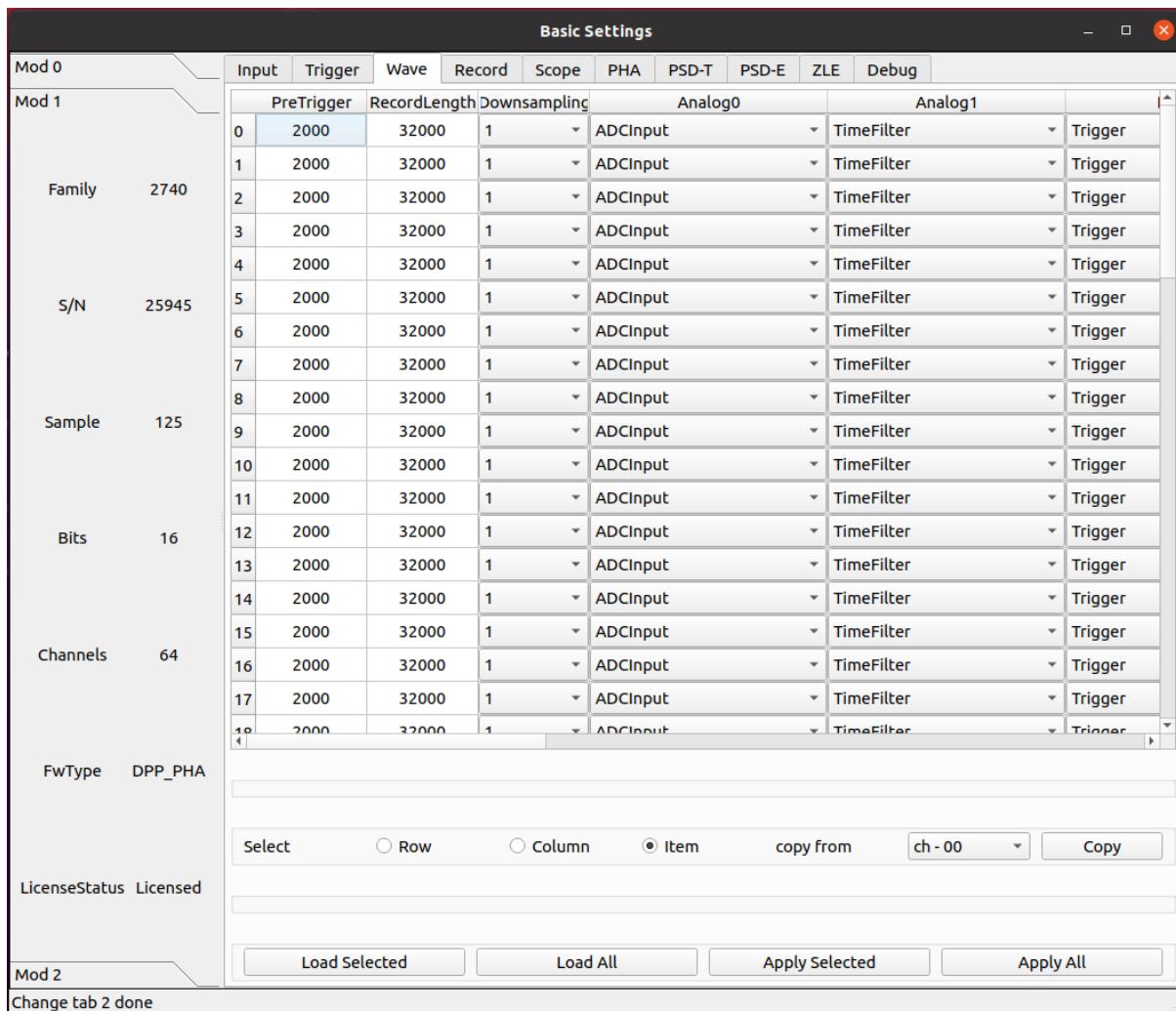
- 1
  - x1
- 2
  - x2
- 4
  - x4
- 8
  - x8

### 参数 Analog0/1


- ADCInput
  - ADC input probe
- TimeFilter
  - Time Filter probe
- EnergyFilter
  - Energy Filter probe
- EnergyFilterBaseline
  - Energy Filter Baseline
- EnergyFilterMinusBaseline
  - [Energy Filter –Baseline] probe

### Digital0/1/2/3

- Trigger
  - Trigger probe
- TimeFilterArmed
  - Time Filter Armed probe
- ReTriggerGuard
  - ReTrigger Guard probe
- EnergyFilterBaselineFreeze
  - Energy Filter Baseline Freeze probe
- EnergyFilterPeaking

- Energy Filter Peaking probe
- **EnergyFilterPeakReady**
  - Energy Filter Peak Ready probe
- **EnergyFilterPileupGuard**
  - Energy Filter Pile Up Guard probe
- **EventPileUp**
  - Event Pile Up probe
- **ADCSaturation**
  - ADC Saturation probe
- **ADCSaturationProtection**
  - ADC Saturation Protection probe
- **PostSaturationEvent**
  - Post Saturation Event probe
- **EnergyFilterSaturation**
  - Energy Filter Saturation probe
- **AcquisitionInhibit**
  - Acquisition Inhibit probe

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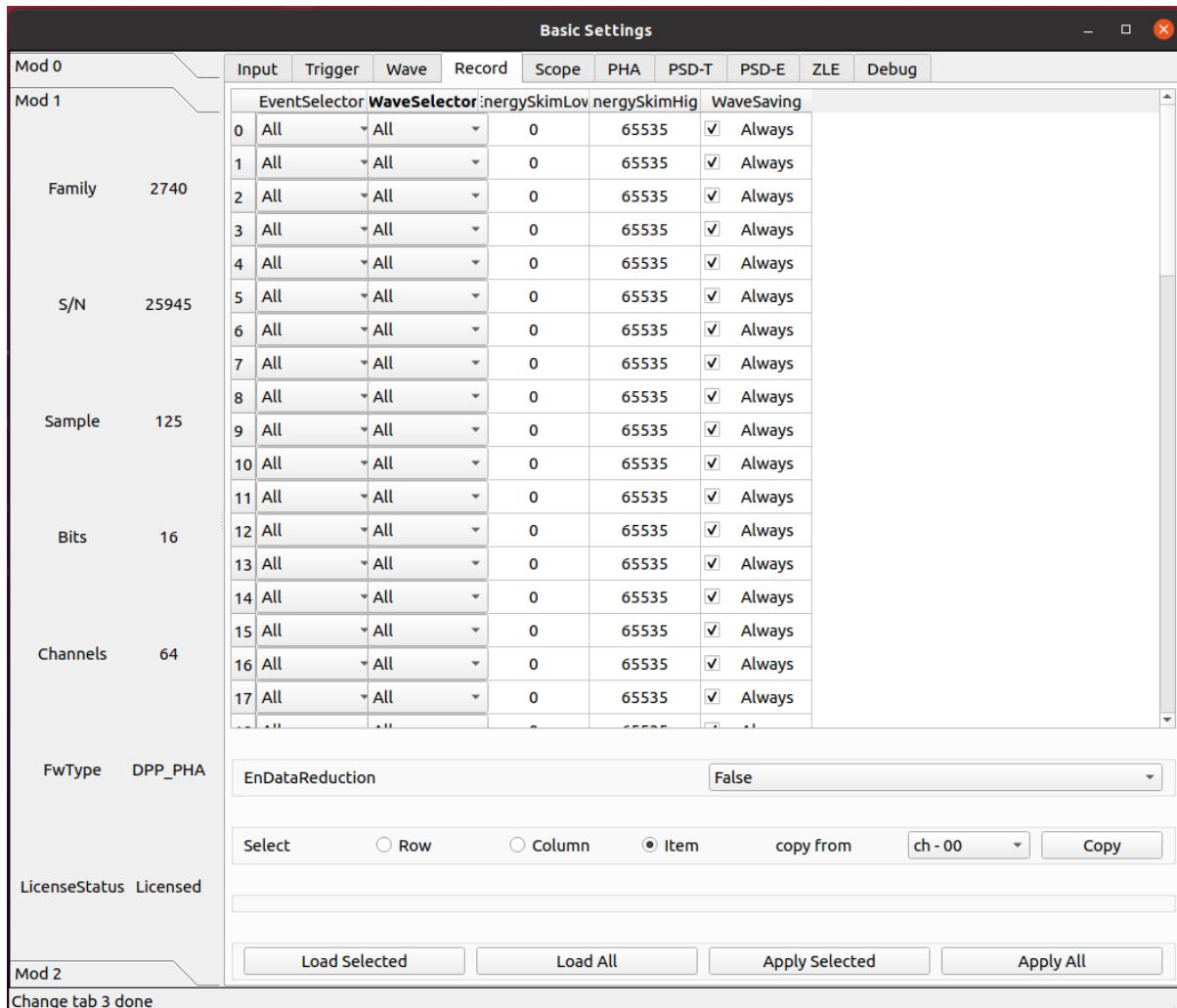### 参数 EventSelector

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- All
  - All events are saved
- PileUp
  - Only pileup events are saved
- EnergySkim
  - Save only the events in the Energy Skim range

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- **Always**
  - Waveforms are always saved
- **OnRequest**
  - Waveforms are saved on request

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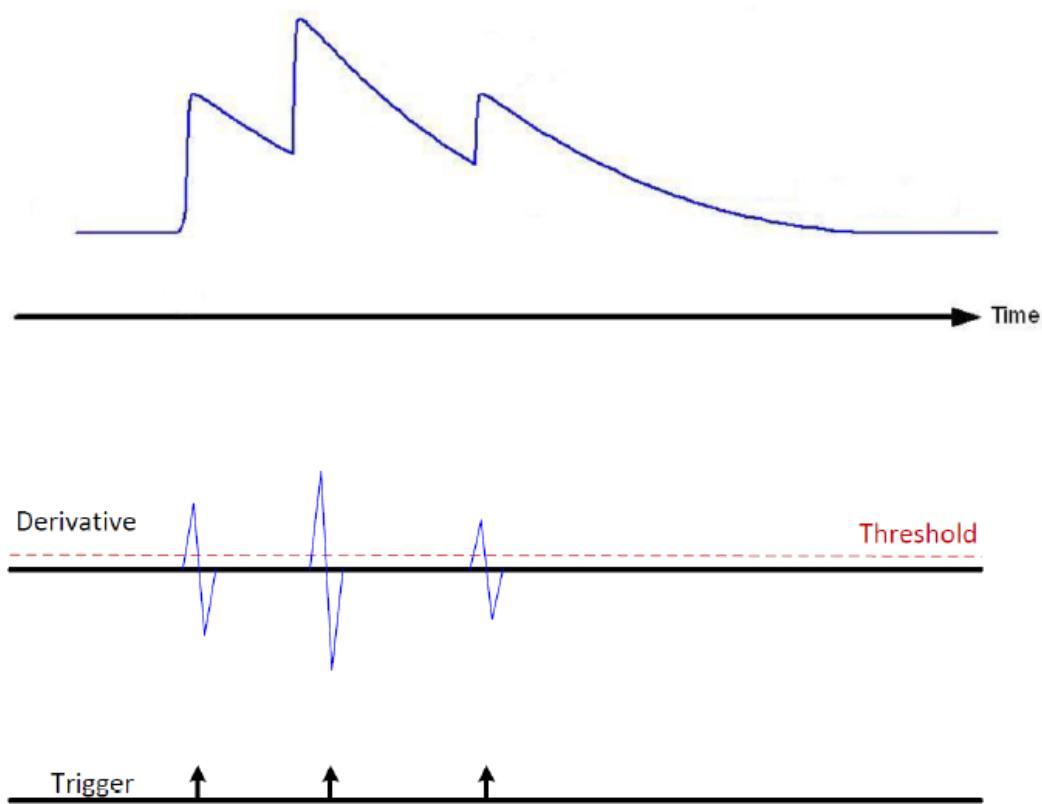
- **Fixed**
  - Baseline fixed at 0
- **VeryLow**
  - Baseline samples for average = 16
- **Low**
  - Baseline samples for average = 64
- **MediumLow**
  - Baseline samples for average = 256
- **Medium**


- Baseline samples for average = 1024
- **MediumHigh**
  - Baseline samples for average = 4096
- **High**
  - Baseline samples for average = 16384

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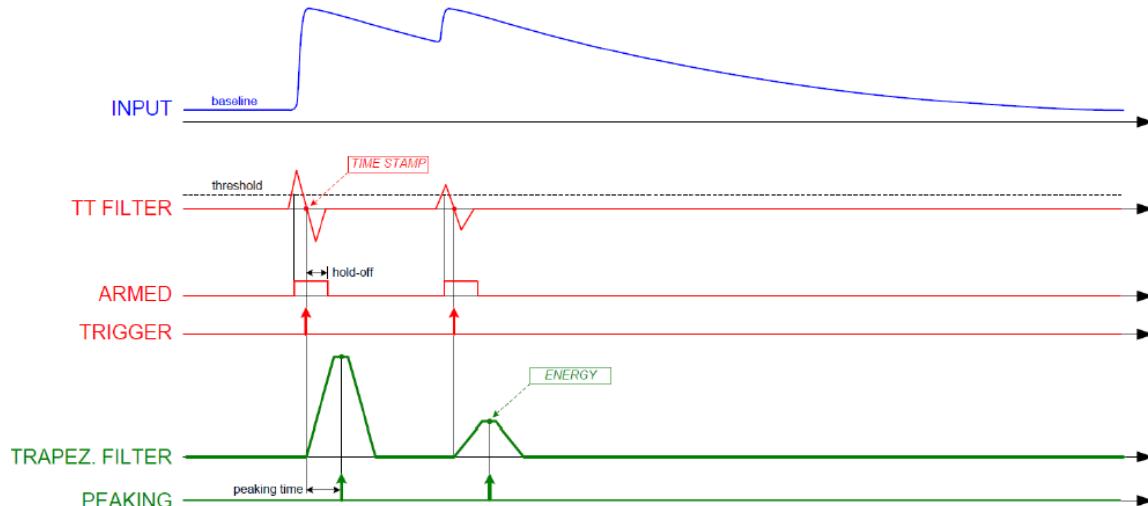- **OneShot**
  - 1 sample
- **LowAVG**
  - 4 samples
- **MediumAVG**
  - 16 samples
- **HighAVG**
  - 64 samples

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
- Long Gate probe
- **PileUpTrigger**
  - Pile Up Trigger probe
- **ShortGate**
  - Short Gate probe
- **ChargeOverRange**
  - Integrated Charge Over Range probe
- **ADCSaturation**
  - ADC Saturation probe
- **ADCInputNegativeOverthreshold**
  - ADC Input Negative Overthreshold probe

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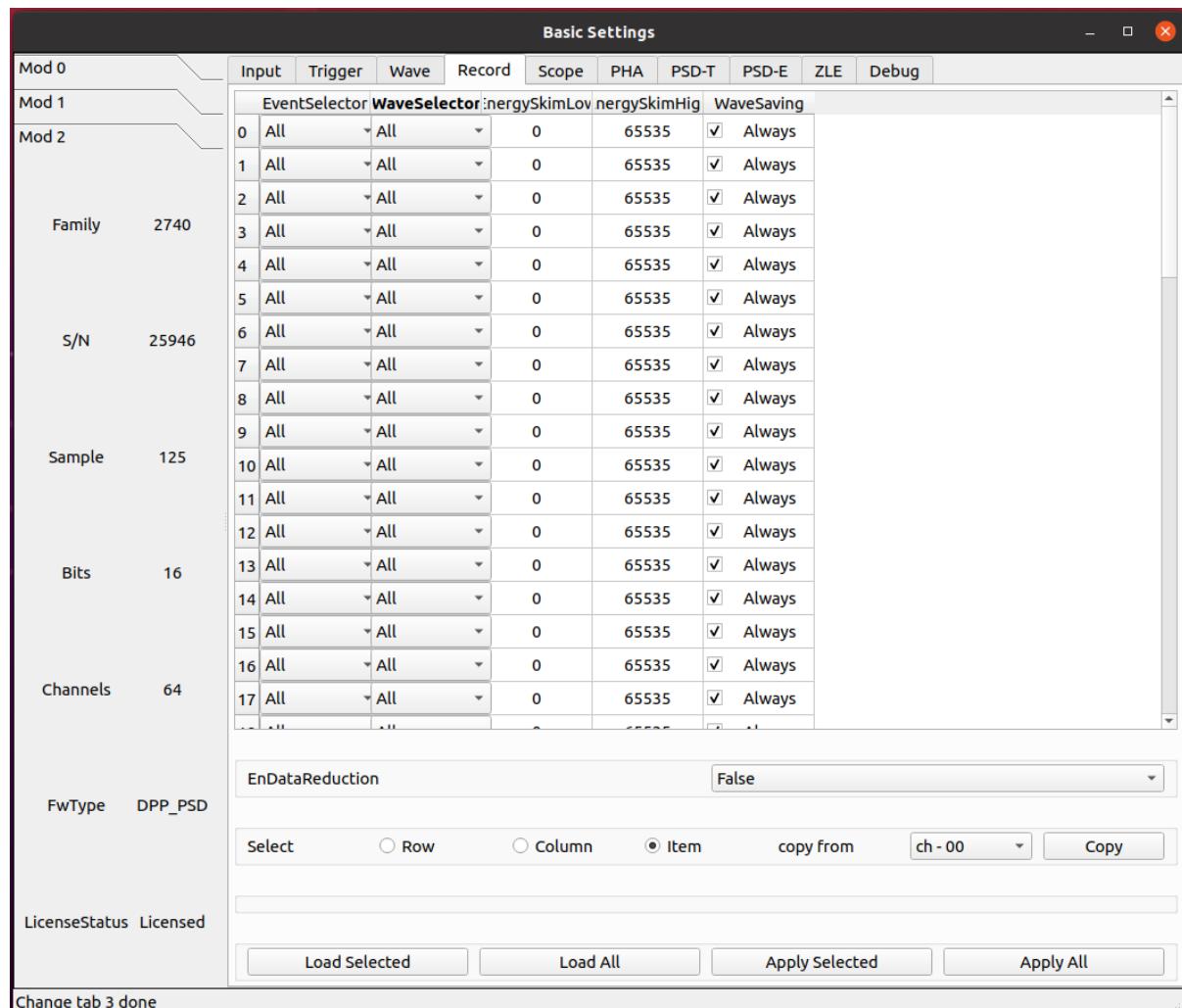### 参数 EventSelector

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- **All**
  - All events are saved
- **PileUp**
  - Only pileup events are saved
- **EnergySkim**
  - Save only the events in the Energy Skim range

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- **All**
  - All waves are saved
- **PileUp**
  - Only pileup waves are saved
- **EnergySkim**
  - Save only waves in the Energy Skim range

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
- 1
  - Smoothing is disabled.
- 2
  - Smoothing is done averaging 2 samples.
- 4
  - Smoothing is done averaging 4 samples.
- 8
  - Smoothing is done averaging 8 samples.

- **16**

- Smoothing is done averaging 16 samples.

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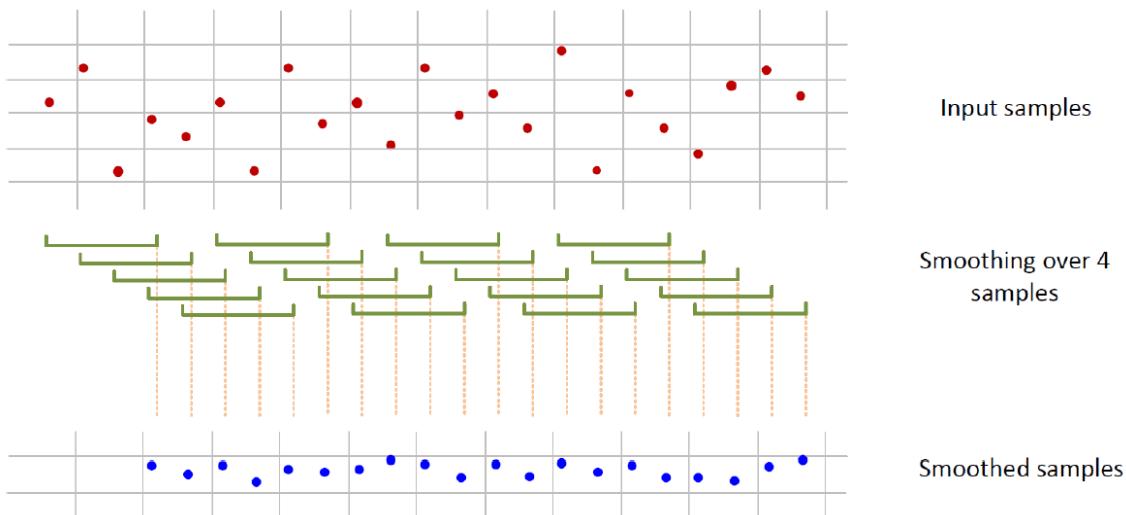- **Enabled**
  - Smoothing factor is enabled for the time filter.
- **Disabled**
  - Smoothing factor is disabled for the time filter.

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
- **LeadingEdge**
  - Set the Leading Edge discriminator
- **CFD**
  - Set the Constant Fraction discriminator

## 参数 CFDDelay

discriminates events based on a CFD signal. The digital CFD signal has been implemented in the classical way except for the input signal inversion. The input waveform is first inverted, then attenuated by a factor  $f$  equal to the desired timing fraction (see CDFraction) of full amplitude, then the signal is inverted again and delayed by a time  $d$  equal to the time it takes the pulse to rise from the constant fraction level to the pulse peak; the latest two signals are summed to produce a bipolar pulse, the CFD, and its zero crossing –corresponding to the fraction  $f$  of the input pulse– is taken as the trigger time.

The delay of the CFD signal can be defined by the user. The TriggerThreshold is then referred to the CFD itself, and the threshold crossing arms the event selection. The trigger fires at the zero crossing of the derivative signal itself.

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
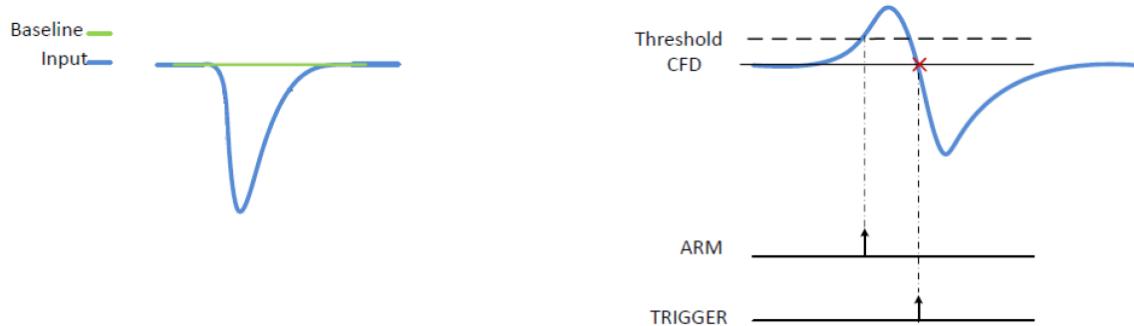
- **Disabled**
  - Trigger hysteresis mechanism is disabled.
- **Enabled**
  - Trigger hysteresis mechanism is enabled.

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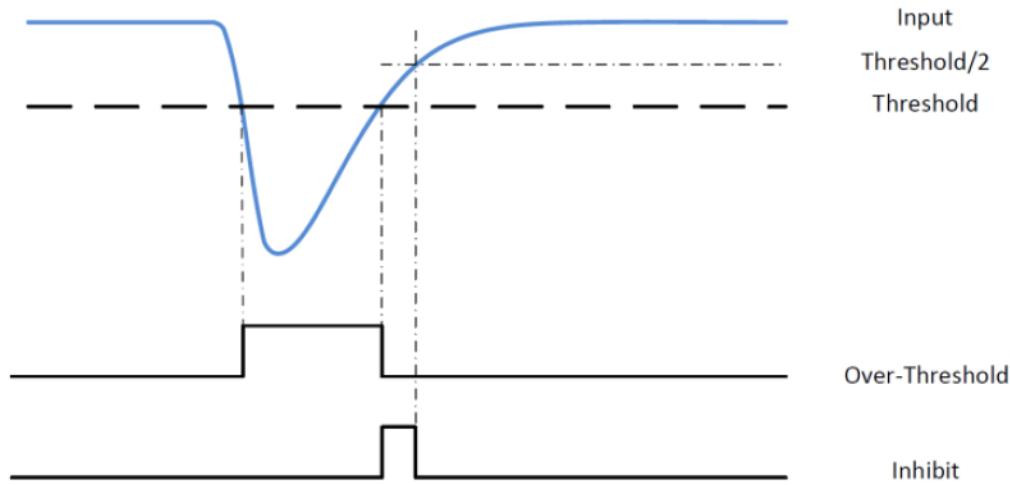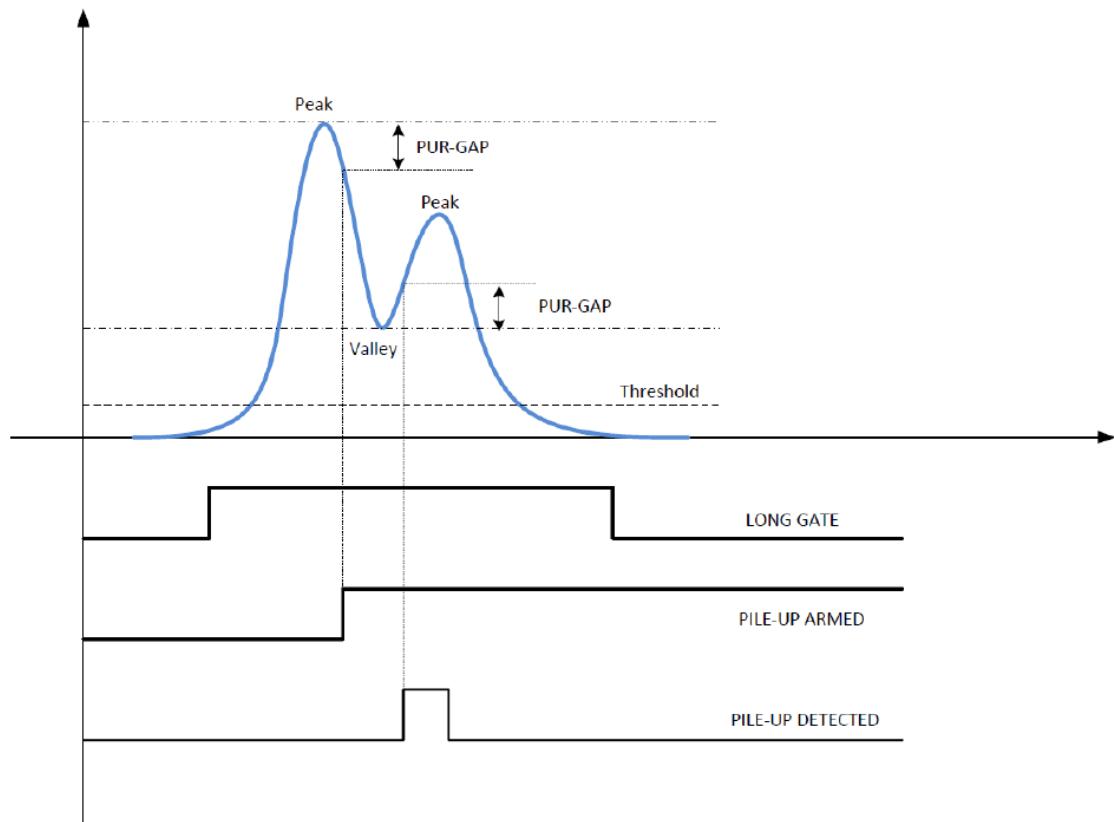
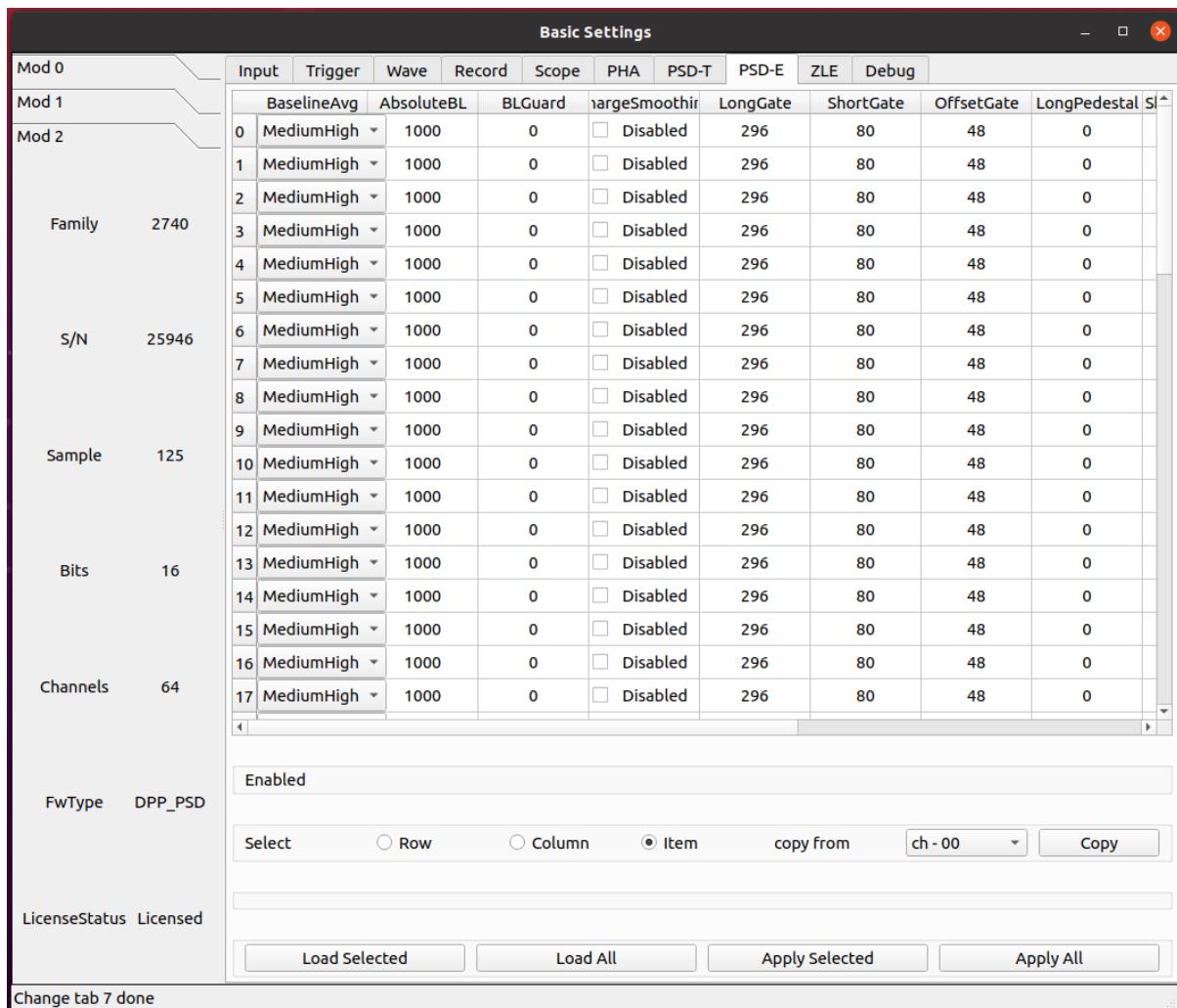
- **Fixed**
  - Baseline fixed at AbsoluteBaseline value
- **Low**
  - Baseline samples for average = 16
- **MediumLow**
  - Baseline samples for average = 64
- **MediumHigh**
  - Baseline samples for average = 256

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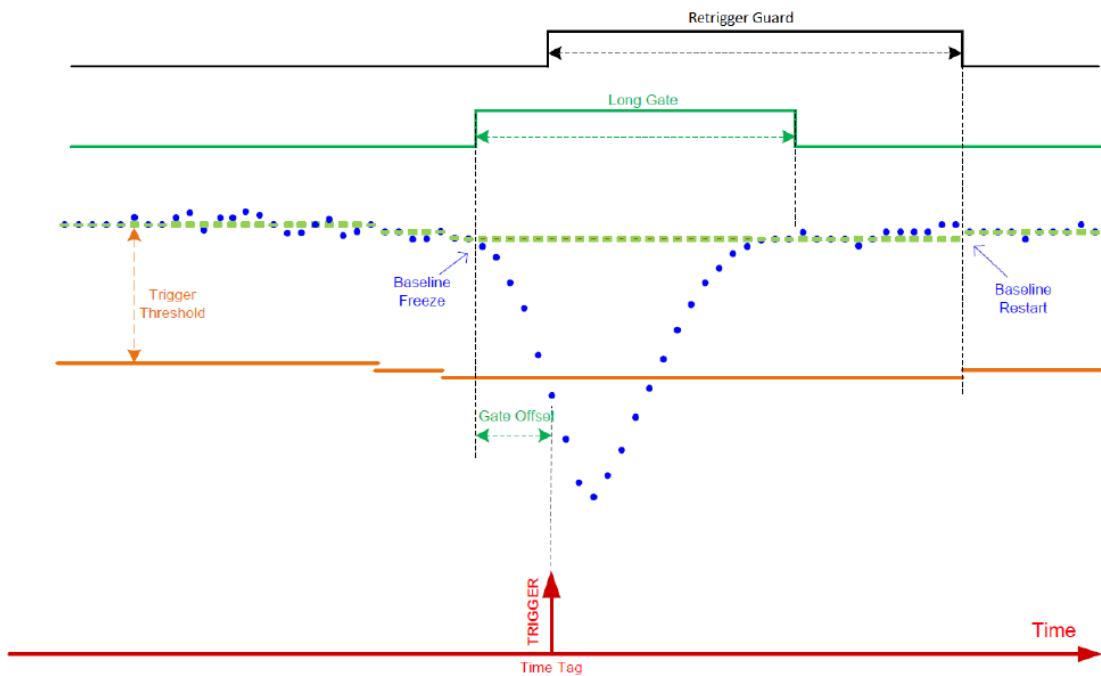- **Enabled**
  - Smoothing factor is enabled in the charge evaluation
- **Disabled**
  - Smoothing factor is disabled in the charge evaluation

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
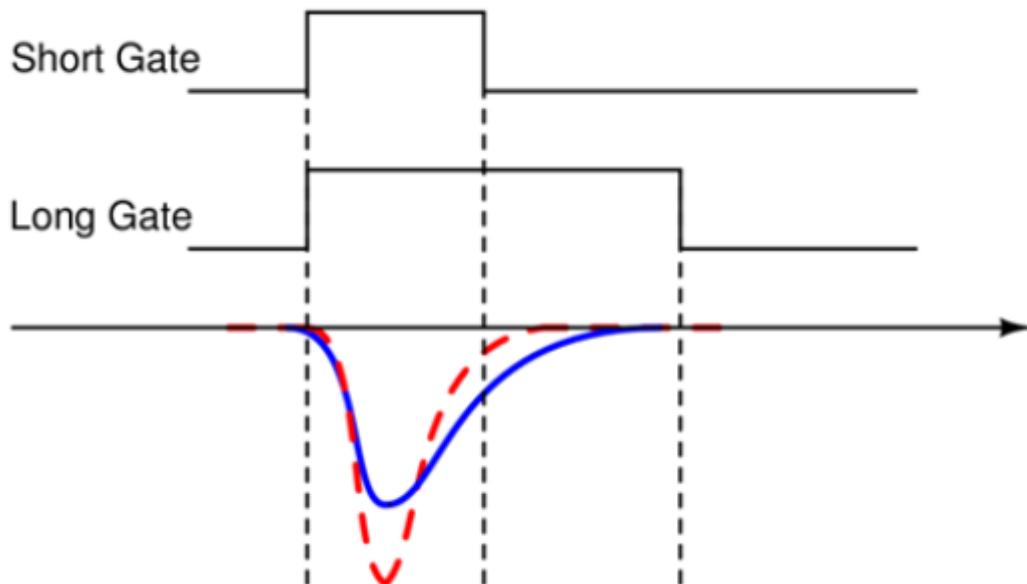
- **x1**
  - Charge value is multiplied x1.
- **x4**
  - Charge value is multiplied x4.
- **x16**
  - Charge value is multiplied x16.
- **x64**
  - Charge value is multiplied x64.
- **x256**
  - Charge value is multiplied x256.

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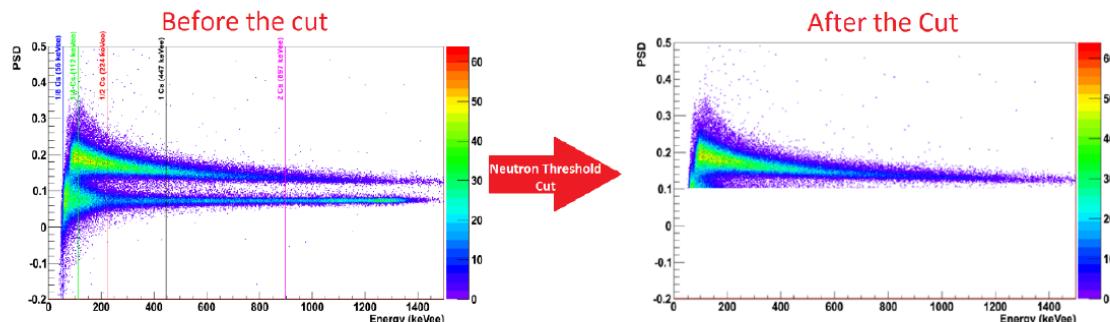- **Disabled**
  - Neutron rejection for events is disabled.
- **Enabled**
  - Neutron rejection for events is enabled.

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
- **Disabled**

- Neutron rejection for waves is disabled.

- **Enabled**

- Neutron rejection for waves is enabled.

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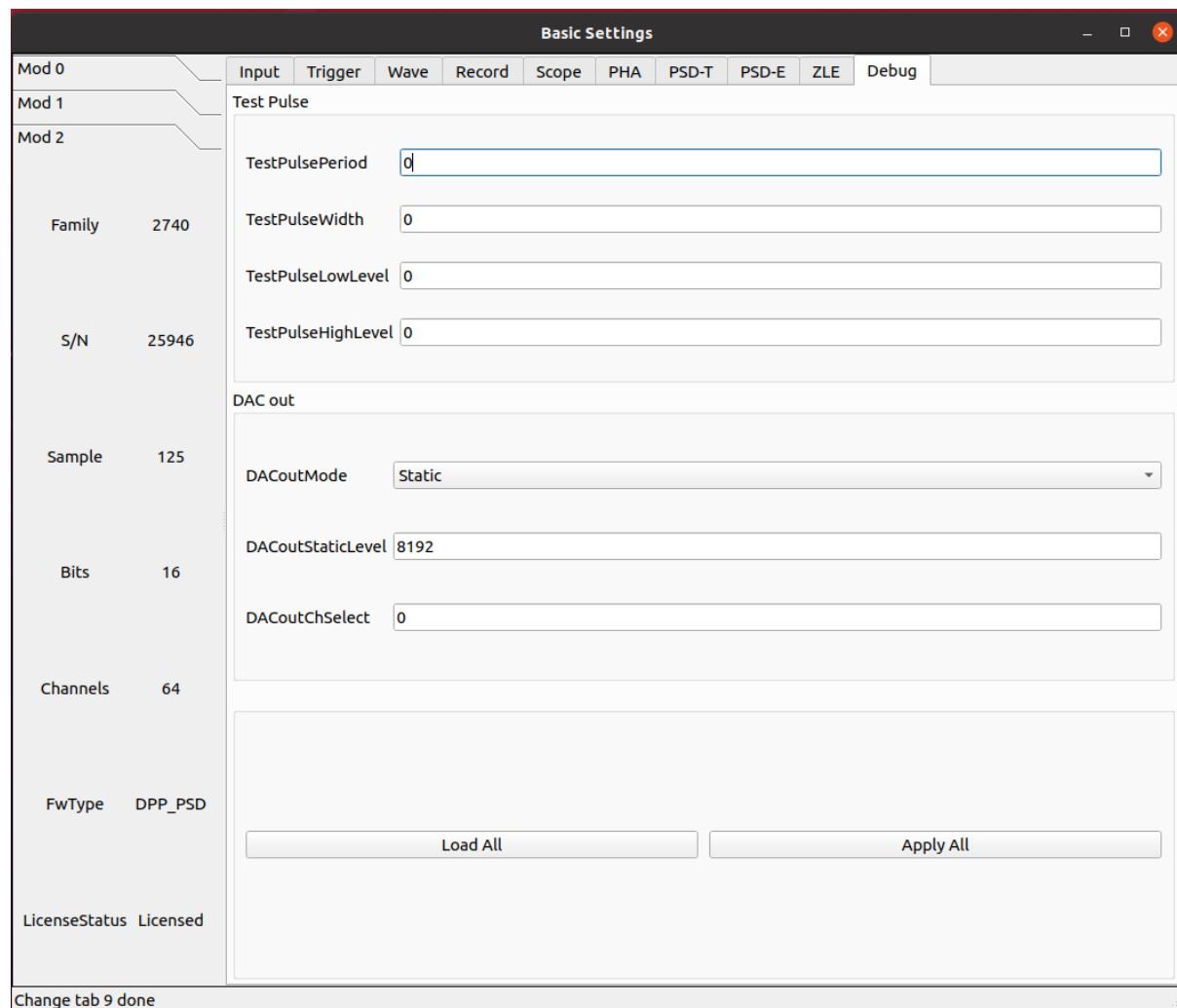- **Static**
  - DAC output stays at a fixed level, given by the DACoutStaticLevel parameter
- **Ramp**
  - The DAC output is driven by a 14-bit counter
- **Sin5MHz**
  - The DAC output is a sine wave at 5 MHz with fixed amplitude
- **Square**
  - Square wave with period and with set by TestPulsePeriod and TestPulseWidth and amplitude between TestPulseLowLevel and TestPulseHighLevel.
- **IPE**
  - Not implemented
- **ChInput**
  - The DAC reproduces the input signal received by one input channel, selected by the DACoutChSelect parameter
- **MemOccupancy**
  - Level of the memory occupancy (not yet implemented)
- **ChSum**
  - The DAC reproduces the “analog” sum of all the digitizer inputs (not yet implemented)
- **OverThrSum**
  - The DAC output is proportional to the number of channels that are currently above the threshold

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### Firmware

This page allows you to upload a new CUP file (firmware image) to your digitizer and select which firmware to run as the default or current version. The 'current' firmware is the one currently running on the digitizer, while the 'default' firmware is the one loaded at startup. You can switch between the loaded firmware on-the-fly, except in cases where changing the communication link is required.

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Model	# Channels	MS/s	# bit	Applications
x2740	64	125	16	64 MCAs for high channel density spectroscopy Good fit for Neutrino and Dark Matter exp. (candidate for Dark Side)
x2751 (*)	16	1000	14	Ultra-fast detectors with ps timing applications
x2725/x2730 (*)	16	250/500	14	Medium-fast detectors Sub-ns timing combined with high energy resolution
x2724 (*)	32	125	16	Spectroscopy & MCA Advanced Front-End (gain, shaping, ...) Semiconductor detector (HPGe, Clover, SDD, ...) Typically connected to charge Sensitive Preamplifier

(\*) Coming soon. Product specifications are subject to change without notice.

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